FitCute®


Kegel Rejuve

User Guide

Contents

English 3
Français 12
Deutsch 21
Nederlands 30
Español 39
Português 48
Italiano 57
Svenska 66
Polski 75
Čeština 84
Русский 93
Українська 102
简体中文 111
繁體中文 120
日本語 129
한국어 138
Manufacturer & EC REP 147

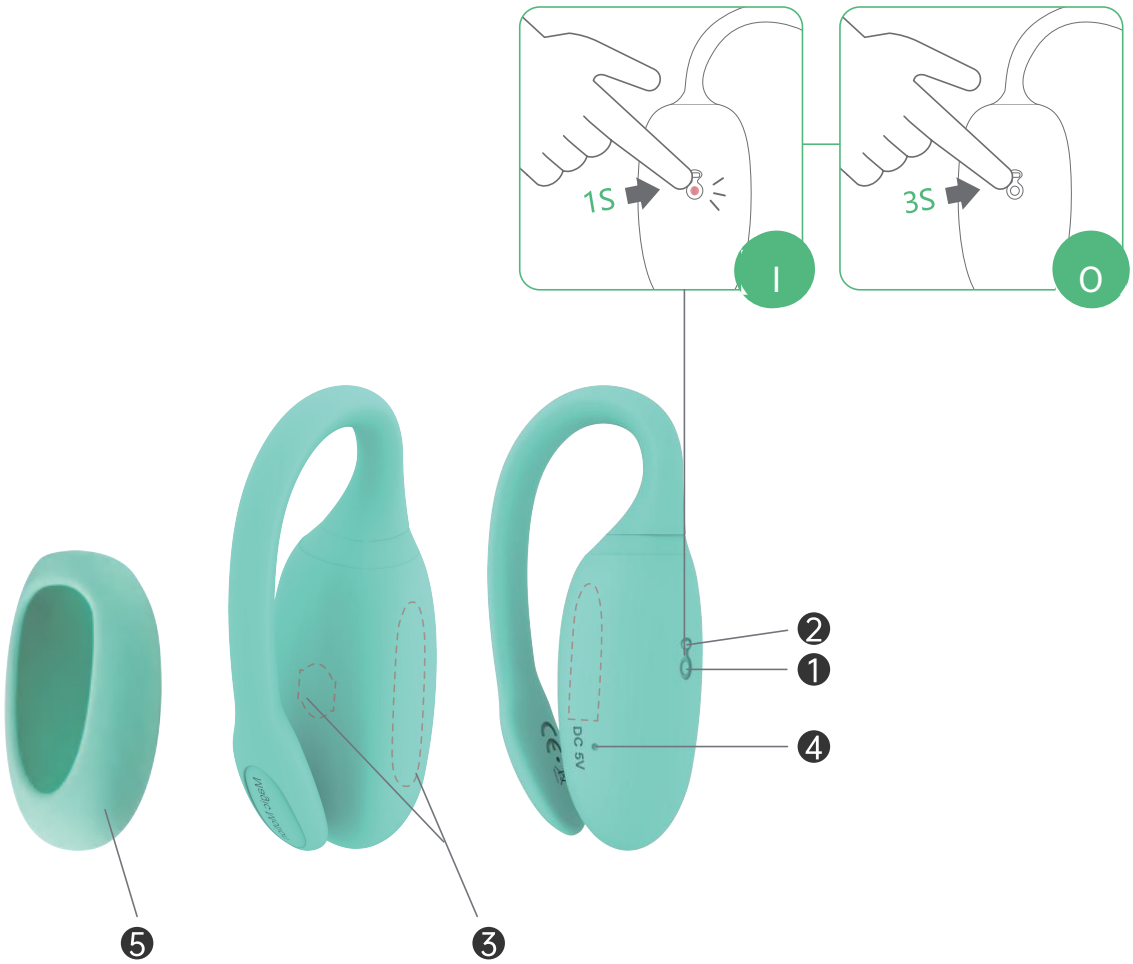1 Power on/off

ON · 1 s

OFF · 3 s


2 LED Light


3 Touch sensor areas

4 Charging port

5 Rejuve Cover

English

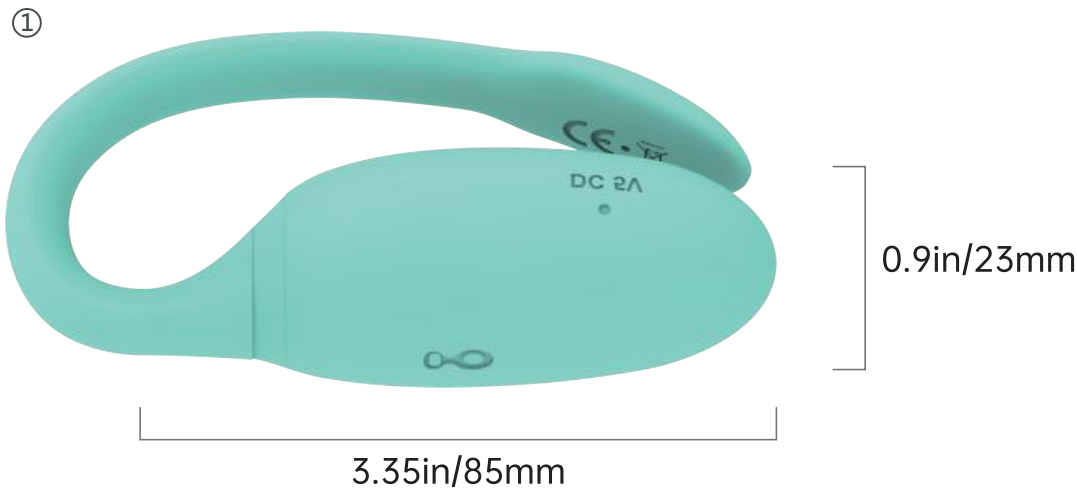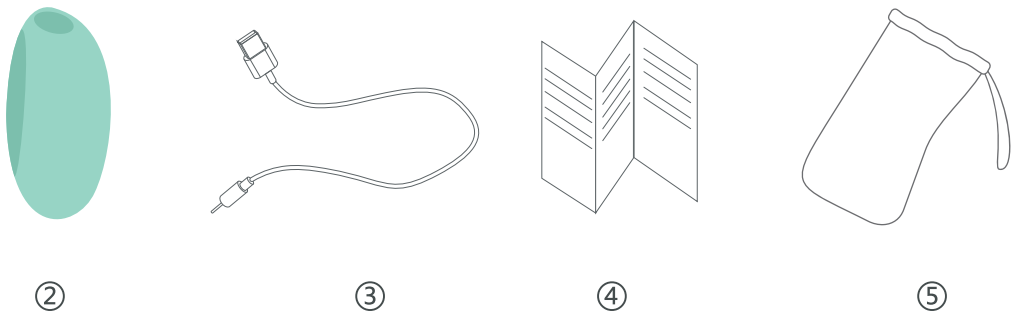1 Kegel Rejuve


2 Rejuve Cover

3 Charging Cable

4 User Guide


5 Storage Bag

Preparation for use

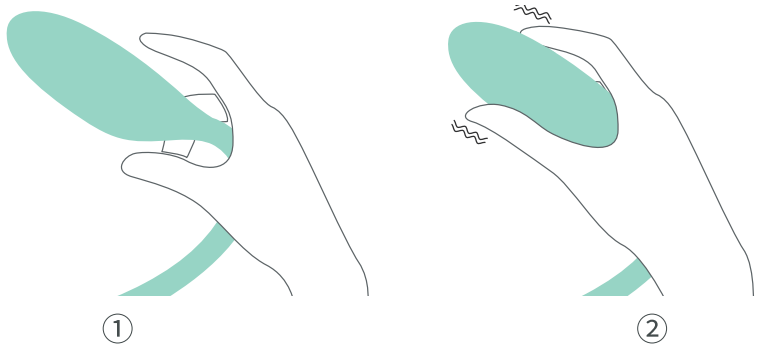Touch the sensor

Find pelvic floor muscle

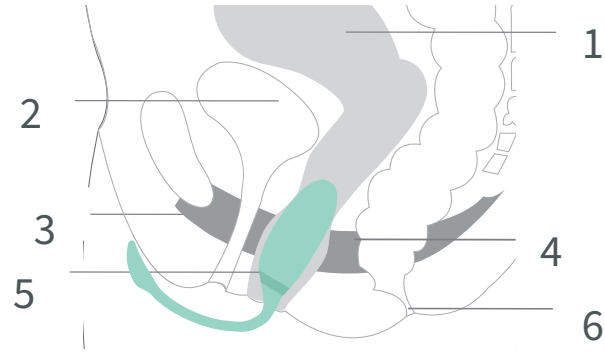1 Uterus

2 Urinary bladder

3 Urethra

4 Pelvic floor muscles


5 Vagina

6 Anus


A water-based lubricant

Using Kegel Rejuve with Rejuve Cover

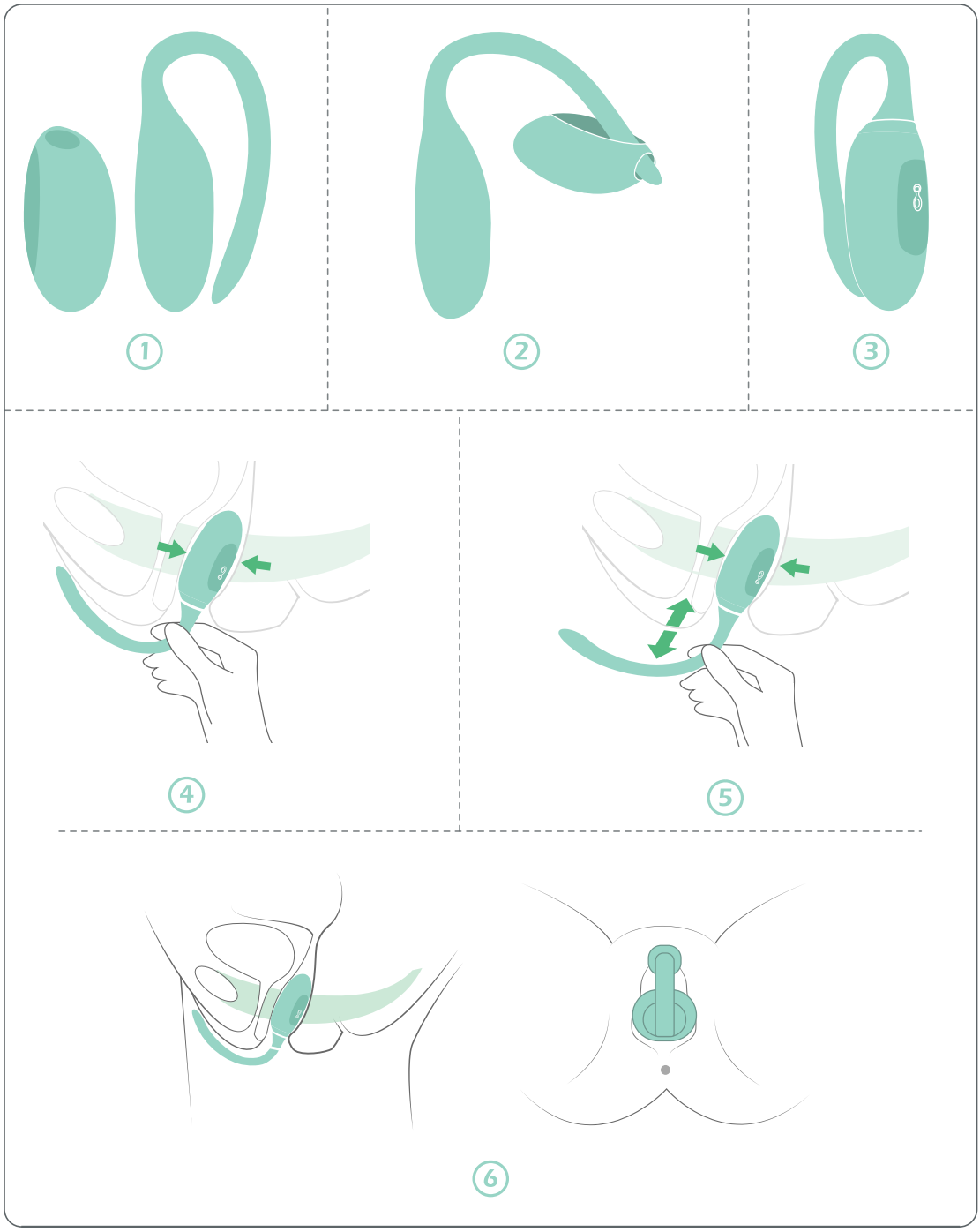Force-Feedback Training Mode

ON · 1 s

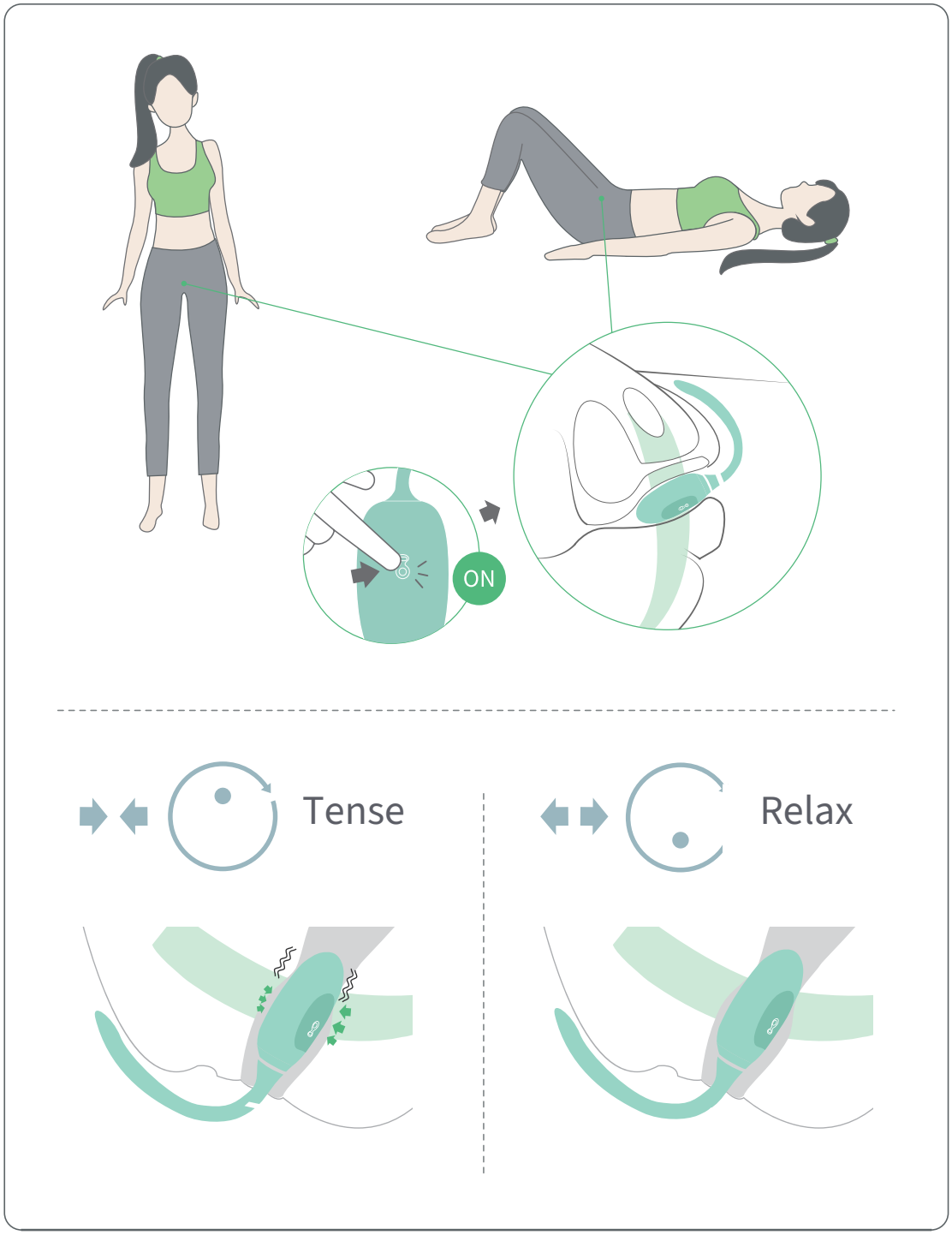Smartphone Workout Mode

Supports iOS & Android smartphones & tablets

1 Download APP

Magic Kegel


Available on the App Store

Google play


2 Connection

Charging

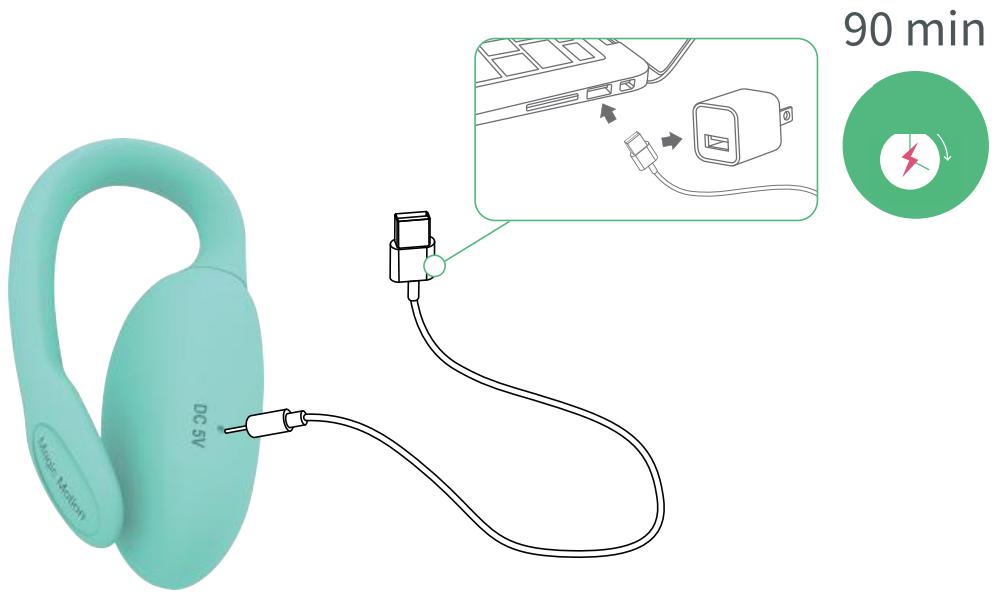Maintenance

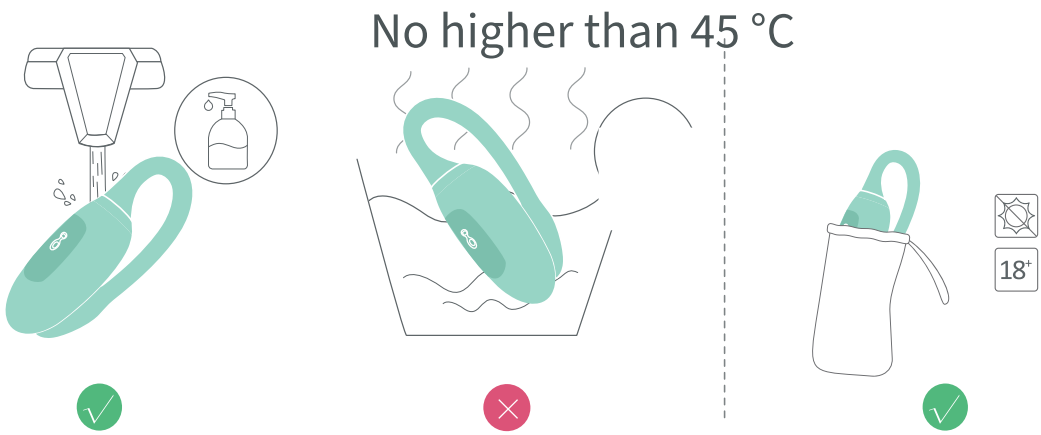18+: adults only

Specification

Model: SMV10

Size: Φ0.91x3.35(6.89) inch / Φ23x85(175)mm

Material: Silicone/ ABS

Charging voltage: DC 5V

Maximum charging current: 300mA

Weight: 27g

Operating time: 150 minutes

Charging time: 90 minutes

Waterproof: IPX4

Operating Frequency Band: 2402–2480 MHz

Maximum RF Transmit Power: 0.95 dBm

Customer Support

For support or warranty questions, email customercare@magicsmotion.com. Include the model and issue; include proof of purchase for warranty questions.


Important Safety Information

Safe Use: For adults aged 18 and over only. Use only as described in this guide. Keep away from children. Risk of injury! Remove piercings or jewellery in the genital area before use. Do not use the toy on irritated or damaged skin. Stop using the toy if you experience pain or feel unwell. Do not use the toy if it is unusually hot, damaged or deformed, or discoloured.

Charging Safety: Use the supplied charging cable and a power supply compatible with the device's rated input. Keep charging contacts dry and the cable and power supply away from water. Do not use during charging. Disconnect when fully charged.

Battery Safety: Do not disassemble, puncture or heat the device or battery, or expose them to fire. Stop use and charging if there is unusual heat, swelling, leakage or damage.

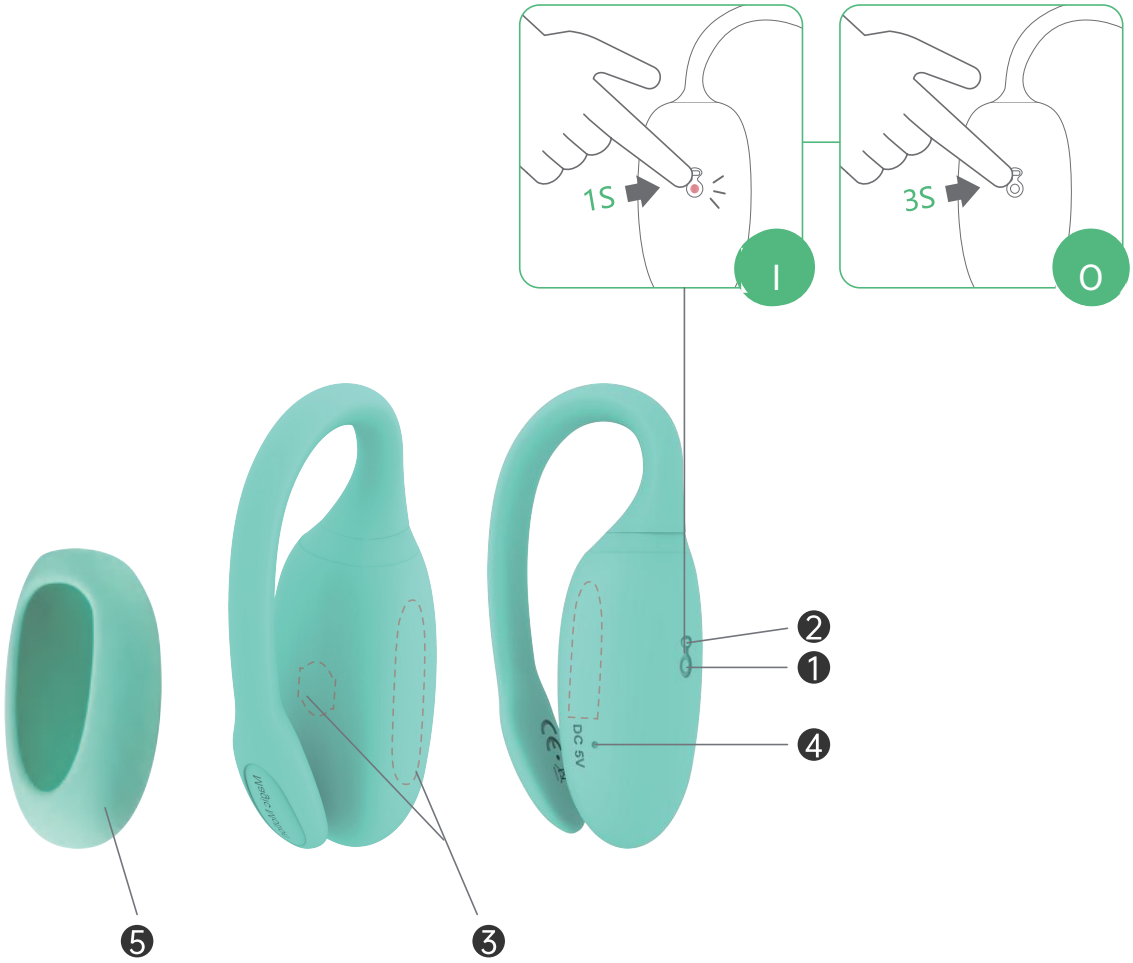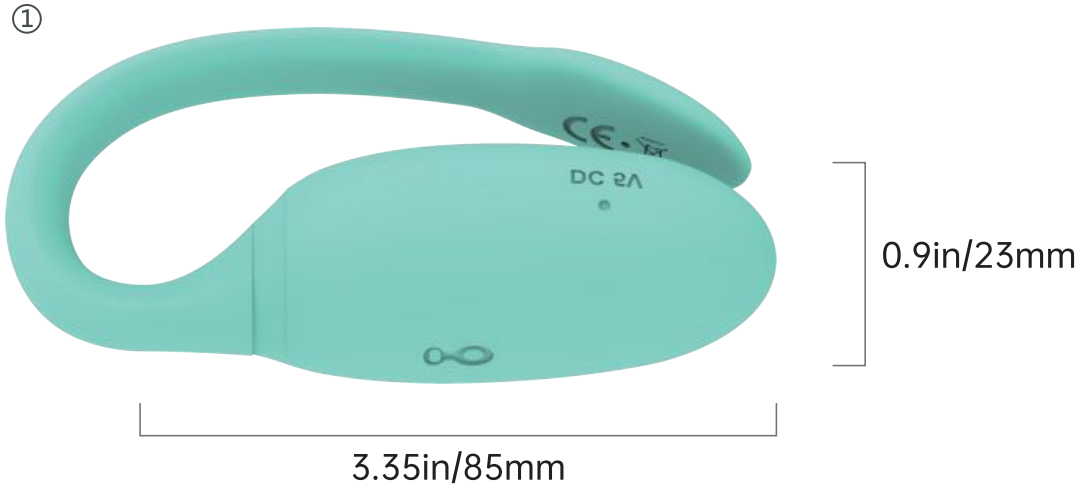
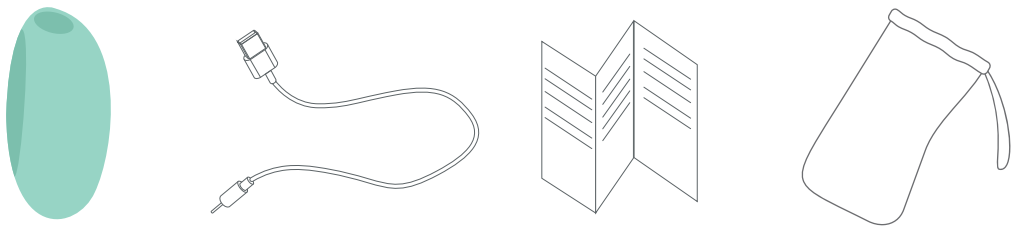
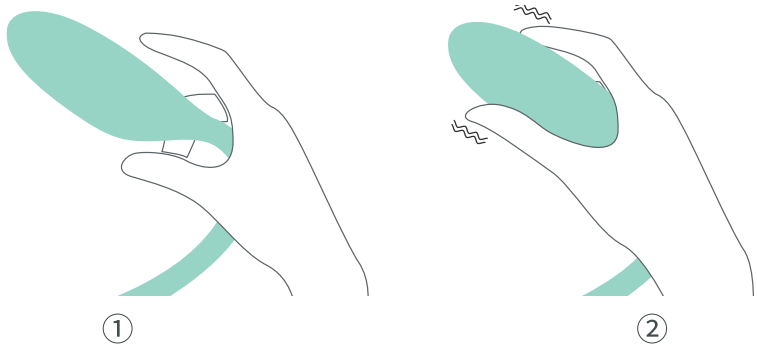
Cleaning, Storage & Disposal: Clean before and after use according to this guide, and dry before storage. Observe the stated IPX4 rating; it does not apply to the charging cable or power supply. Store dry, away from sunlight and children. Take the device and its battery to an appropriate electronics/battery collection point; do not dispose of them with household waste.

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Vtrump Tech (Shanghai) Co., Ltd.