
Magic Motion

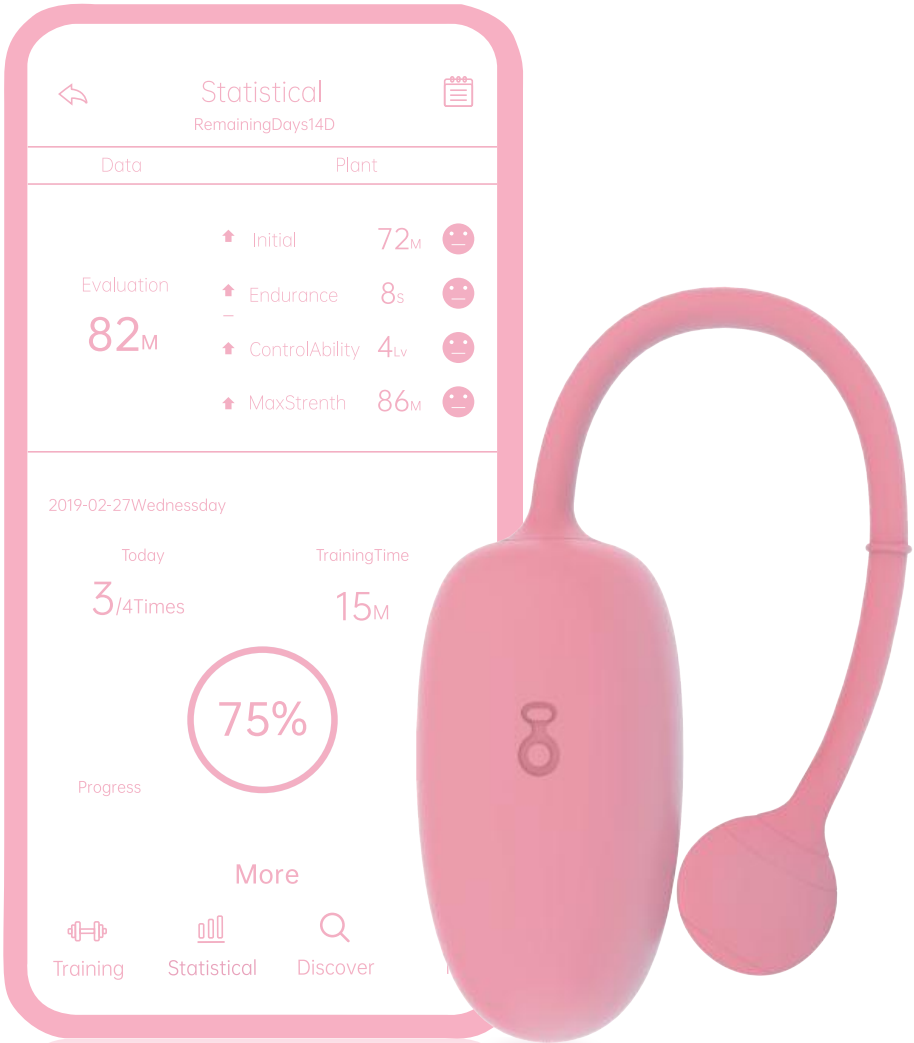
Kegel Coach

User Guide

Contents

English 3
Français 13
Deutsch 23
Nederlands 33
Español 43
Português 53
Italiano 63
Svenska 73
Polski 83
Čeština 93
Русский 103
Українська 113
简体中文 123
繁體中文 133
日本語 143
한국어 153
Manufacturer & EC REP 163

English

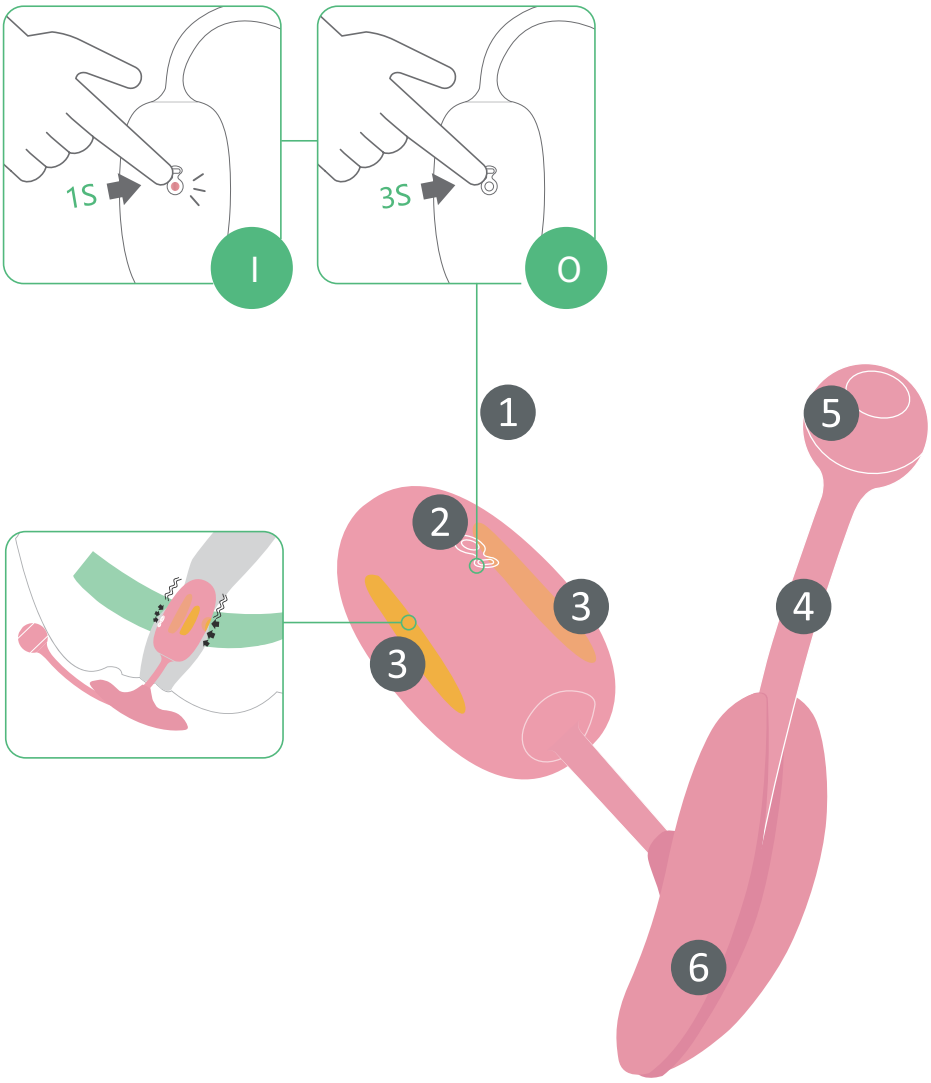ON · 1 s

OFF · 3 s

1 Power on/off

2 LED Light


3 Touch sensor areas


4 Antenna


5 Charging port


6 Anchor

English

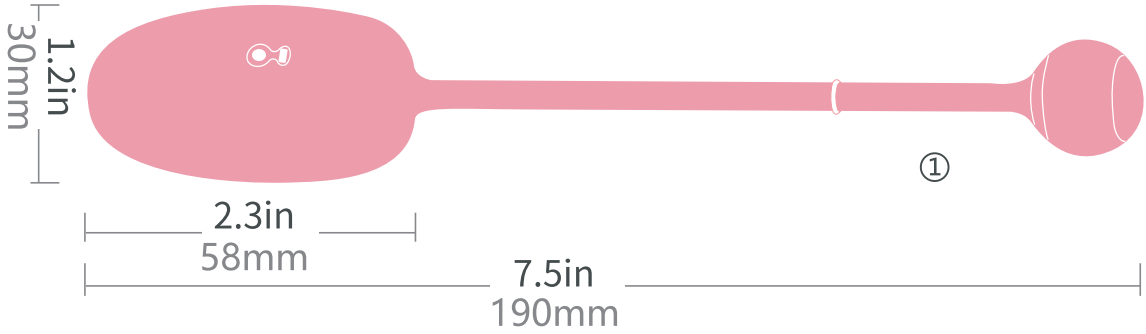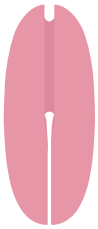②

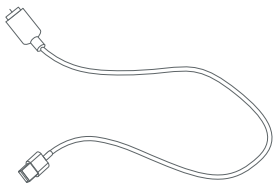③

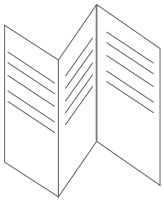④

⑤

1 Kegel Coach


2 Anchor

3 Magnetic USB cable

4 User Guide


5 Storage Bag

Preparation for use

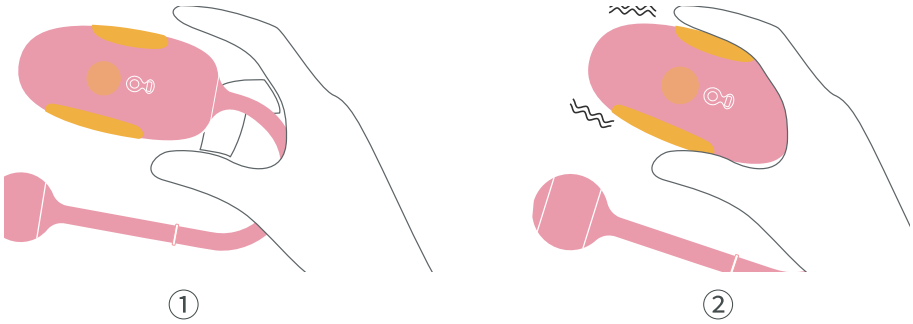Touch the sensor

Find pelvic floor muscle

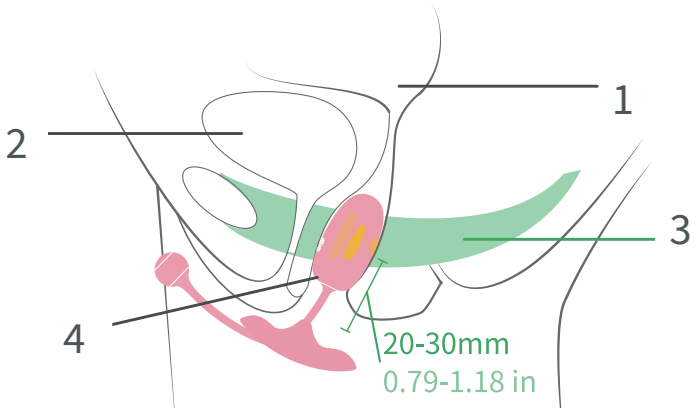
1 Uterus

3 Pelvic floor muscles


2 Urinary bladder

4 Vagina

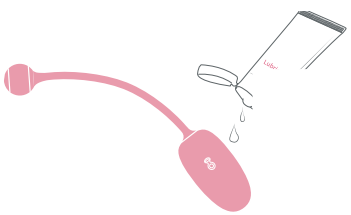
A water-based lubricant

Using Kegel Coach with the anchor

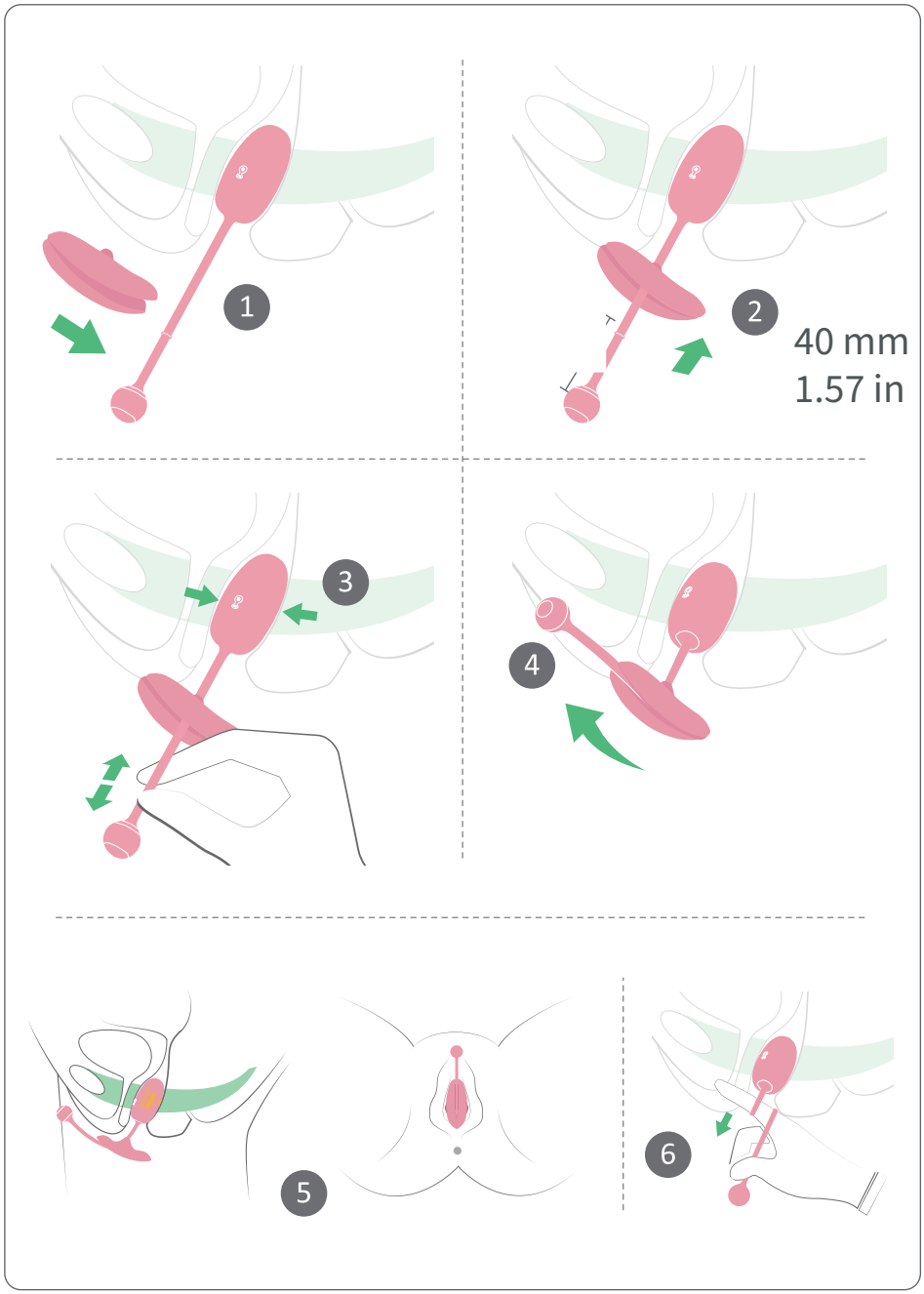Force-Feedback Training Mode

ON · 1 s

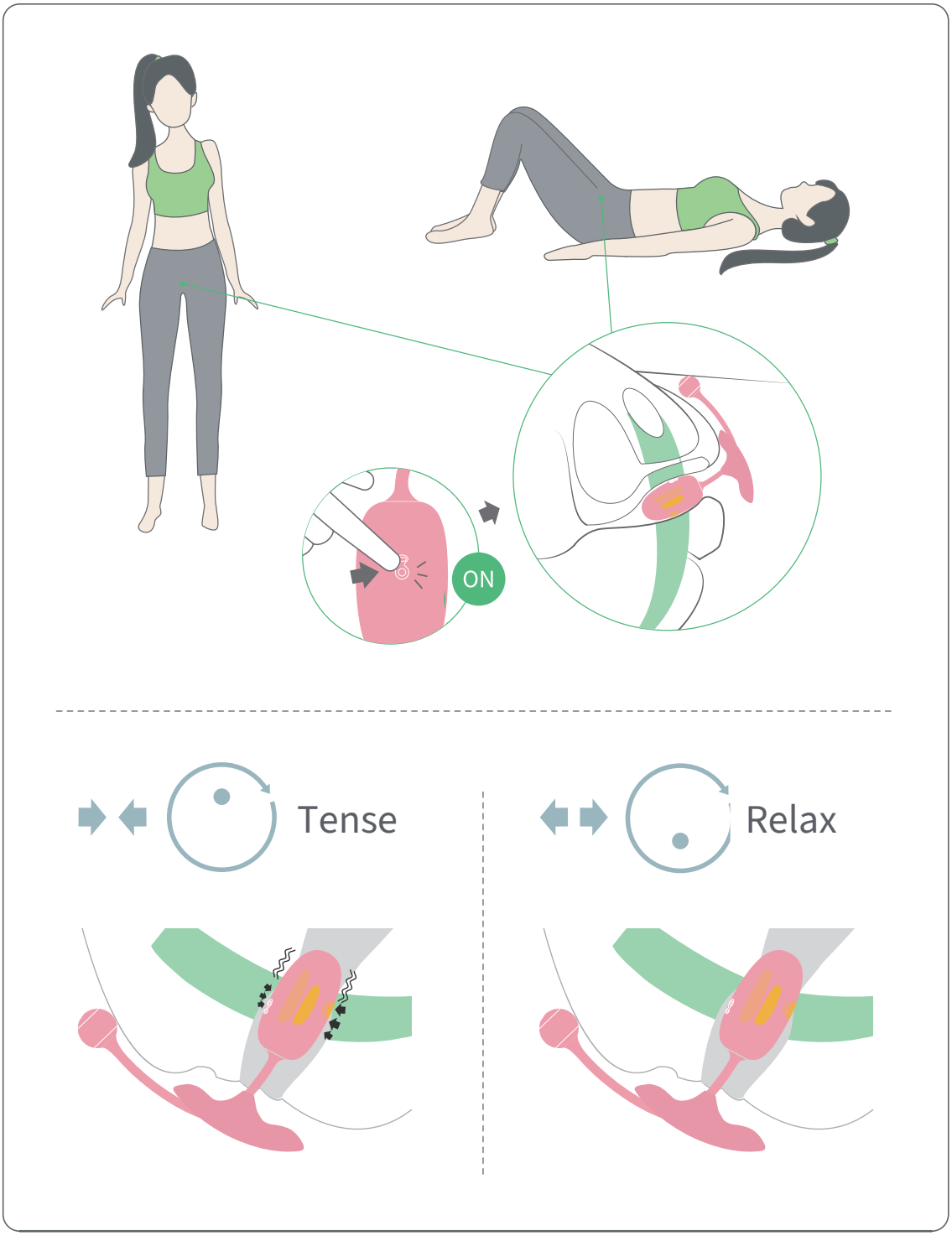Smartphone Workout Mode

Supports iOS & Android smartphones & tablets

1 Download APP


Magic Kegel


Available on the App Store

Google play


2 Connection

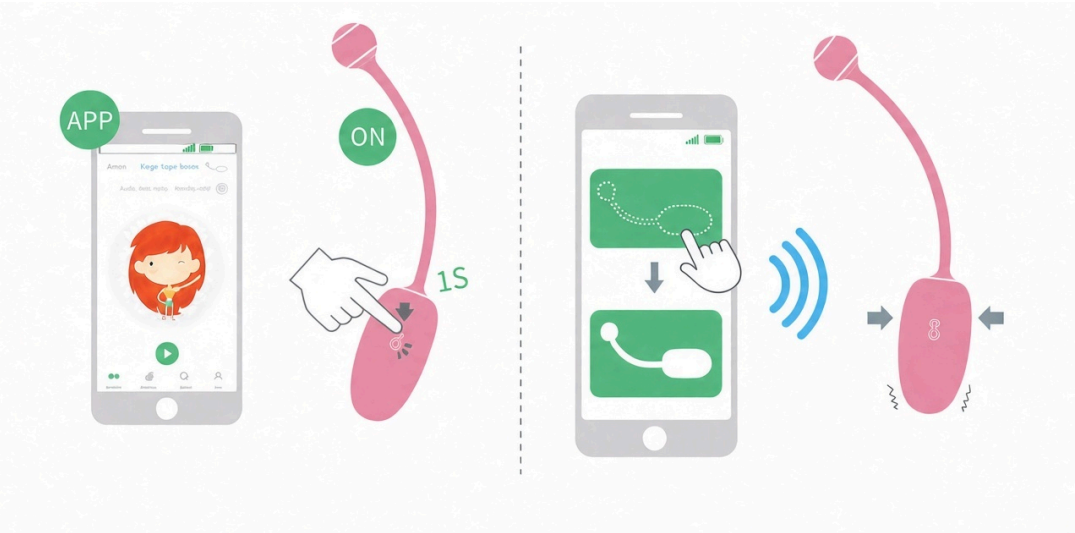Charging

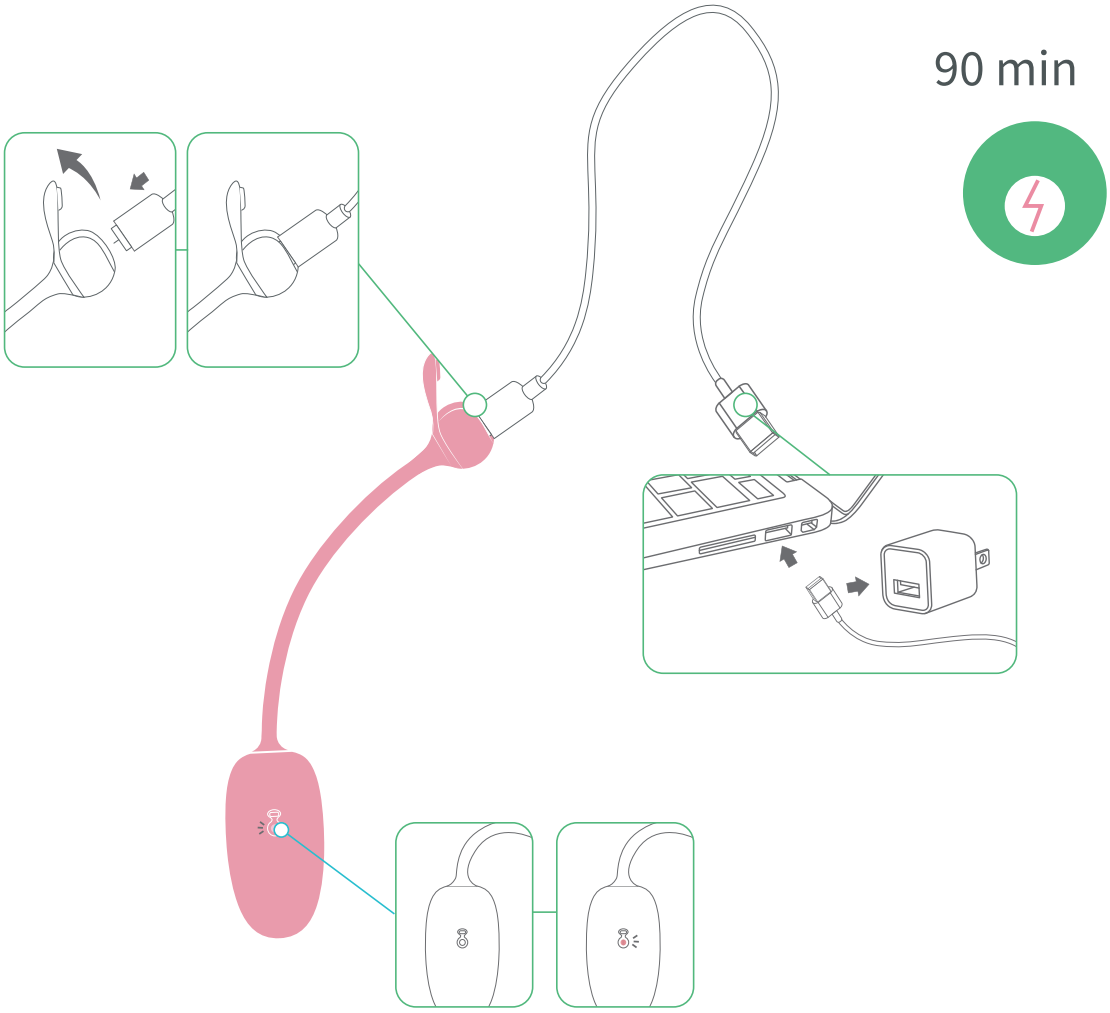Maintenance

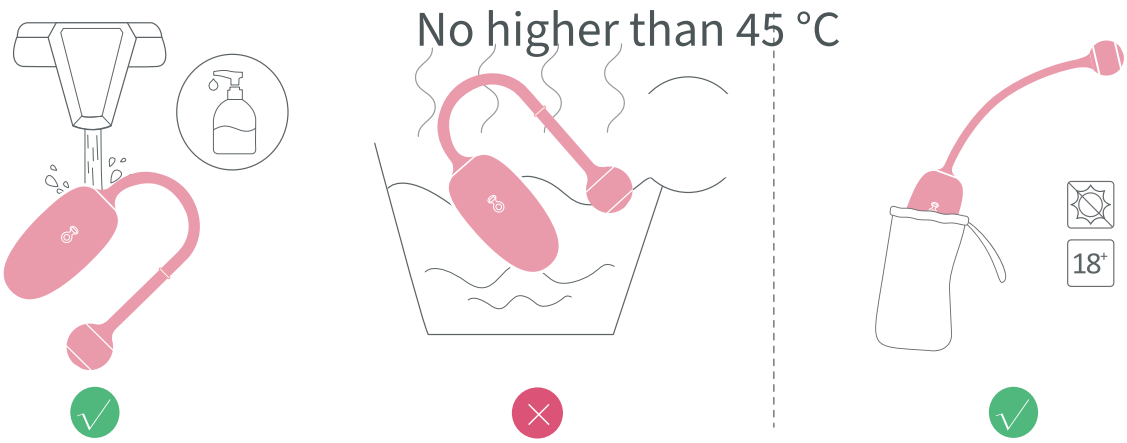18+: adults only

Specification

Model: SMV21

Size: Φ1.2x2.3(7.5) inch

Material: Silicone/ ABS

Charging voltage: DC 5V

Maximum charging current: 300mA

Weight: 36g

Operating time: 150 minutes

Charging time: 90 minutes

Waterproof: IPX6

Operating Frequency Band: 2402–2480 MHz

Maximum RF Transmit Power: 4.61 dBm

Customer Support

For support or warranty questions, email customercare@magicsmotion.com. Include the model and issue; include proof of purchase for warranty questions.

Important Safety Information

Safe Use: For adults aged 18 and over only. Use only as described in this guide. Keep away from children. Risk of injury! Remove piercings or jewellery in the genital area before use. Do not use the toy on irritated or damaged skin. Stop using the toy if you experience pain or feel unwell. Do not use the toy if it is unusually hot, damaged or deformed, or discoloured.

Charging Safety: Use the supplied charging cable and a power supply compatible with the device's rated input. Keep charging contacts dry and the cable and power supply away from water. Do not use during charging. Disconnect when fully charged.

Battery Safety: Do not disassemble, puncture or heat the device or battery, or expose them to fire. Stop use and charging if there is unusual heat, swelling, leakage or damage.

Cleaning, Storage & Disposal: Clean before and after use according to this guide, and dry before storage. Observe the stated IPX6 rating; it does not apply to the charging cable or power supply. Store dry, away from sunlight and children. Take the device and its battery to an appropriate electronics/battery collection point; do not dispose of them with household waste.

Precautions

Do not use device if:

• You have given birth within the last six weeks


• You experience persistent pain when inserting or using device, even with water-based lubricant

• You have a medical condition affecting your pelvic health (for example, abnormal vaginal discharge or severe prolapse)

You can use Rebelle during pregnancy as long as you are not suffering from any complications. Rebelle is an exercise tracker for core muscle strengthening. It must not be used to treat or prevent any medical condition. If you are unsure whether or not you should use device, please consult a qualified medical professional. For hygiene reasons, do not share your device. Do not leave your device in for more than one hour at a time.

Disclaimer

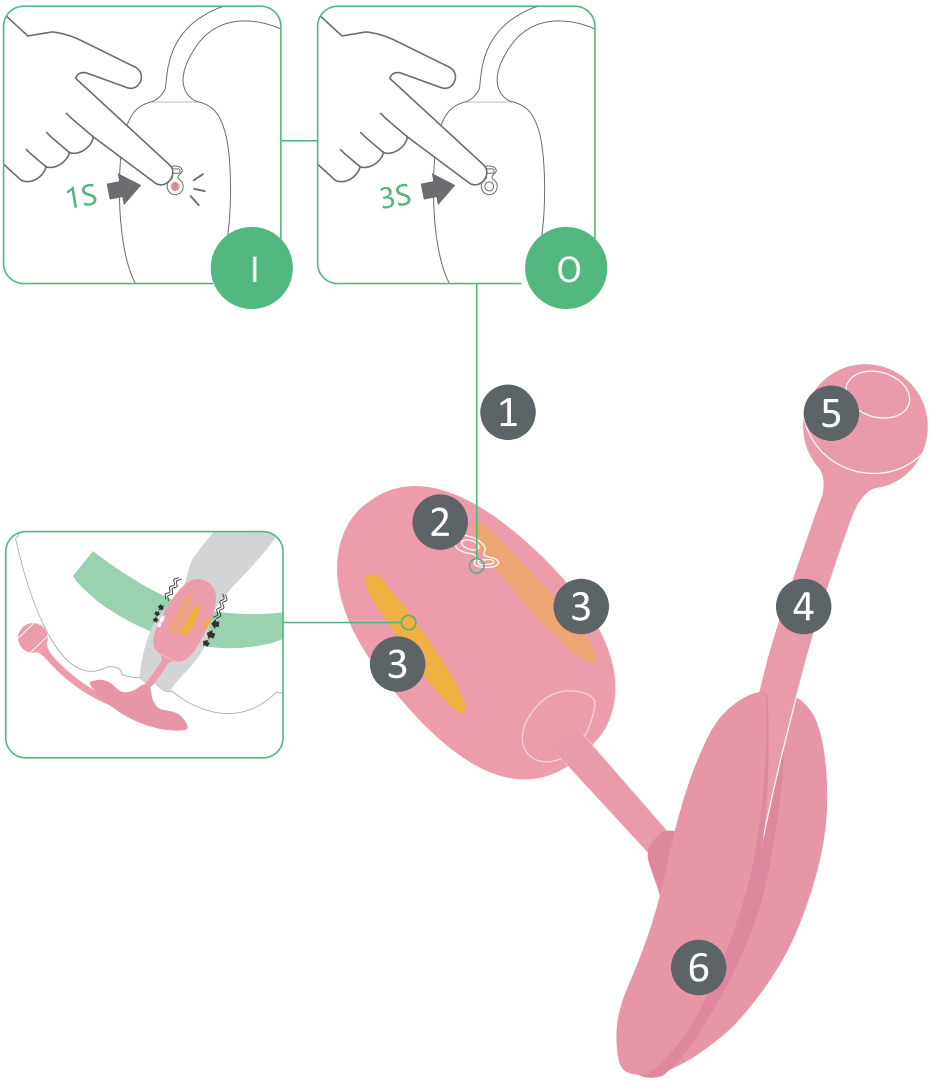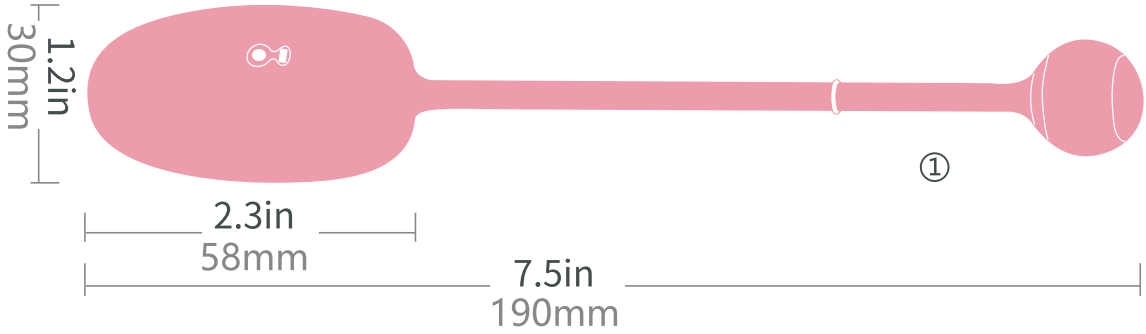
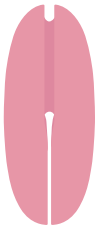
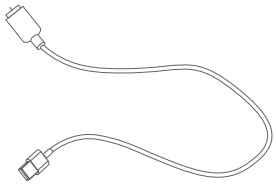
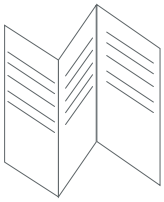
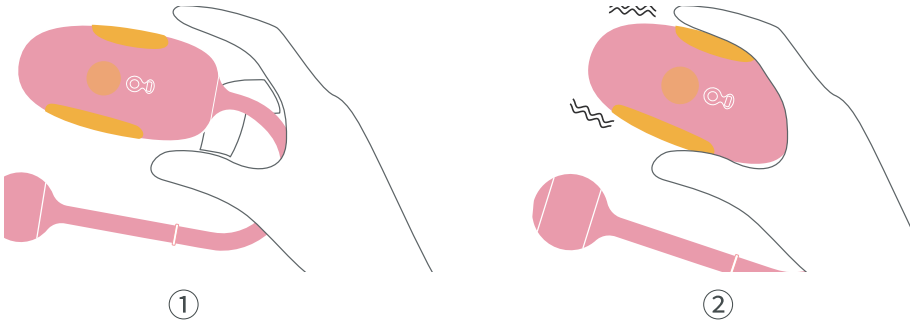
Users of device do so at their own risk. Neither Magic Motion nor its retailers assume any responsibility or liability associated with or for the use of this product. Magic Motion reserves the right to revise and/or update this manual without being obliged to notify any users.

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Vtrump Tech (Shanghai) Co., Ltd.