
Kegel Master

Gen 2

Your Personal Coach Kegel Exercise

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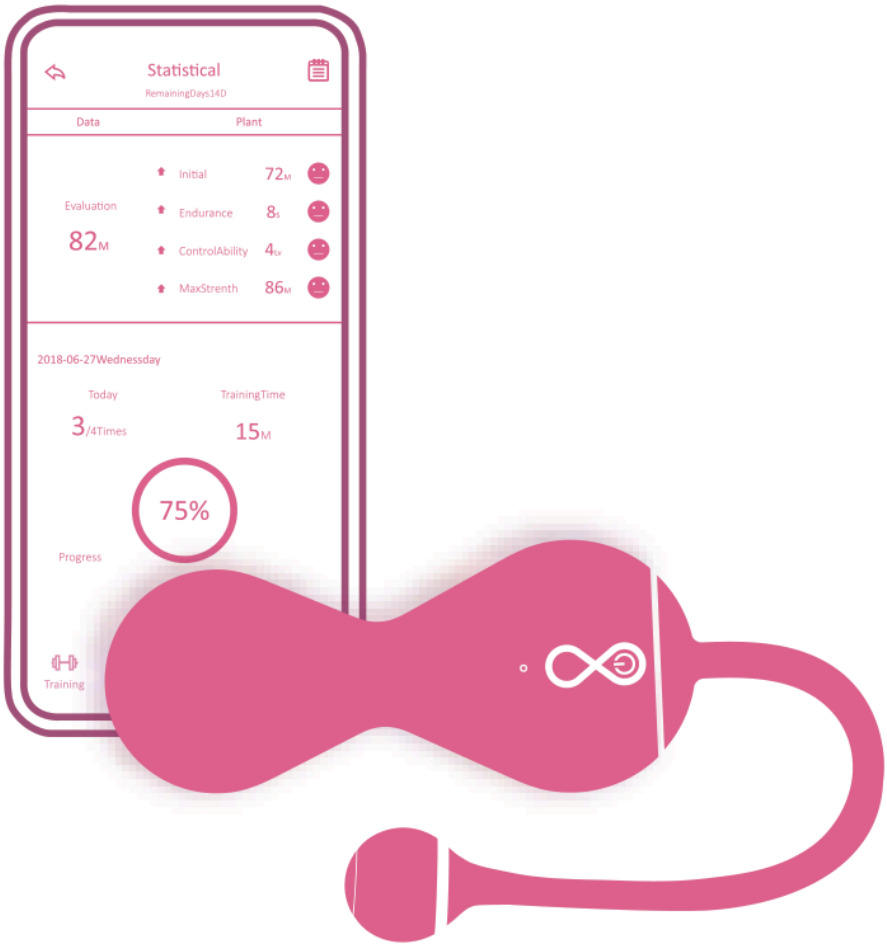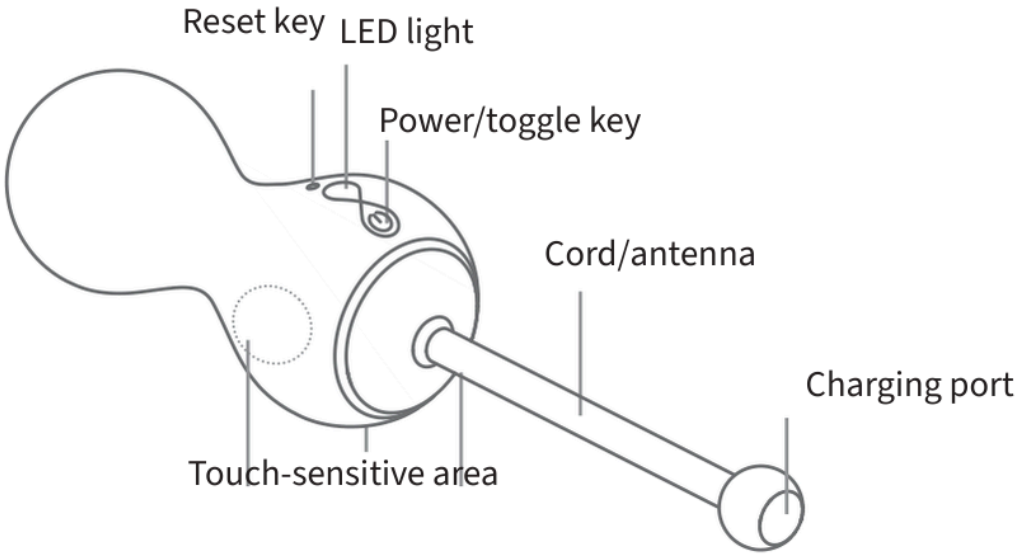
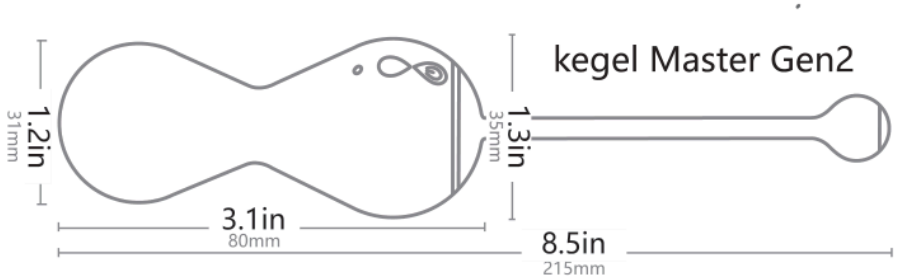
Smartphone Workout Mode

Step1: APP Download

1. Search for "**Magic Kegel**" in IOS App Store or Google Play


2. Scan below QR code to download directly.

3. Download from


https://www.magicsmotion.com/kegel.html

Step2: Run "Magic Kegel" APP


Turn on your
Kegel Master Gen 2 device

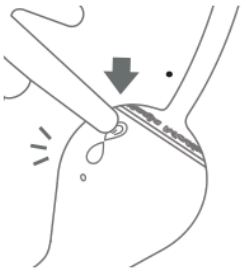Follow the app guide to connect.

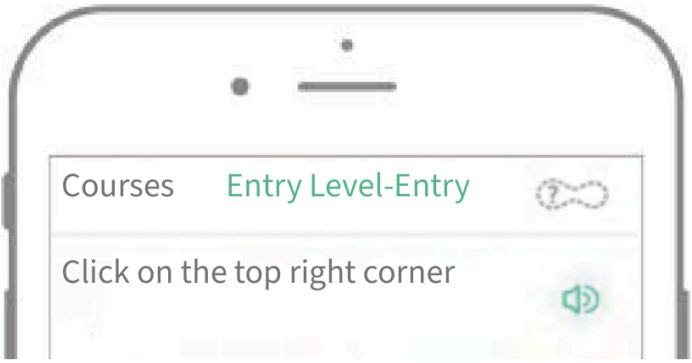The device is connected when this message appears.

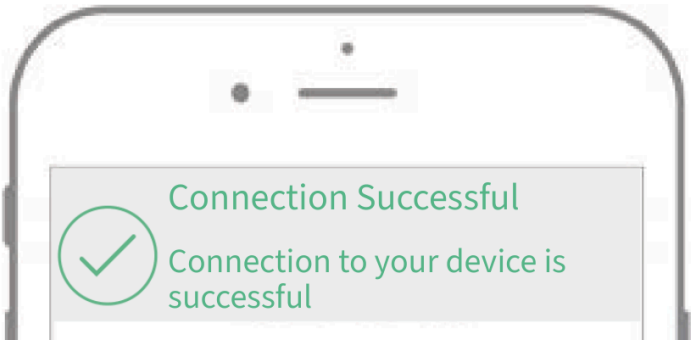Step3: Enjoy Kegel Training


Find a comfortable position for inserting. Many women find lying on the back or standing as most comfortable. Kegel Master Gen 2 LED facing upward, please check below picture. For pushing forward easily, apply a bit water-based lubricant on the top of device or on the condom.

Tips: Please leave the cord of Kegel Master Gen 2 outside the body, at least 4 cm length, to keep good signal without barrier.

Step4: Power Off


The device will turn off after pressing Power Key for 3 seconds, the LED light will turn off to indicate turn-off status.

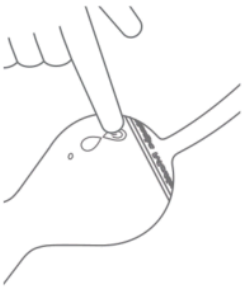Tips: Kegel Master Gen 2 will turn off automatically when it disconnected with App in 5 minutes.


Independent Mode (Without SmartPhone) **- Squeeze for vibration**

Turn on Kegel Master Gen 2 and squeeze pelvic floor muscles to control the strength of vibration of device.

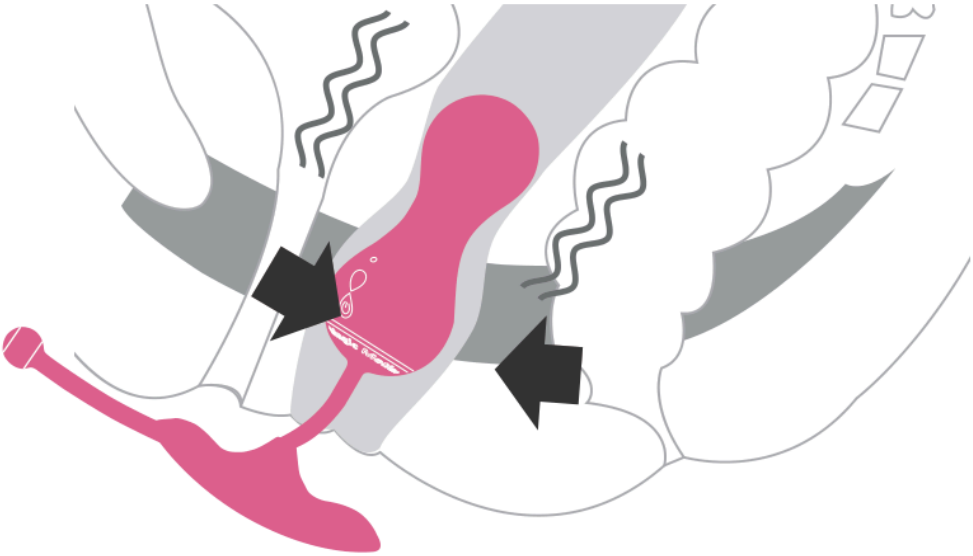
Charging Battery


Step 1 Open the cap of the magnetic charging port for charging.

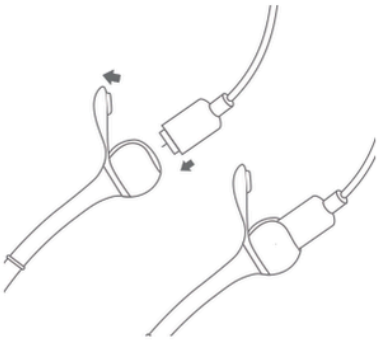Step 2 Connect USB port to smartphone power adapter, laptop, Computer or power charging bank.

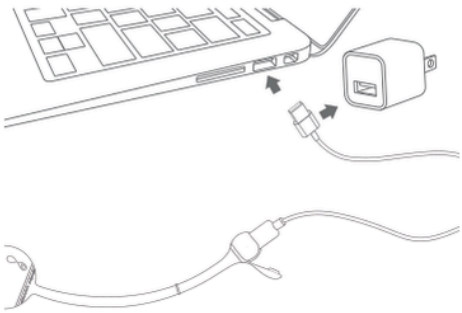Step3 LED light on to show charging normally.The LED light will be off when the charging is finished.

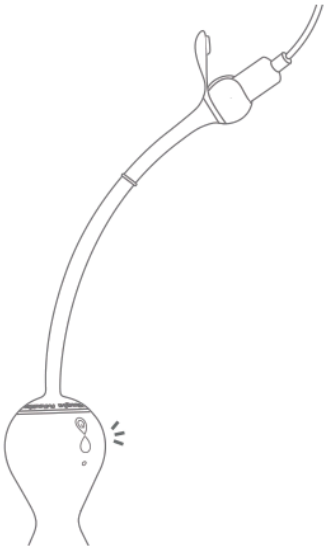Tips: For first use, please charge the battery for 1 hour after the power is fully used.

Cleaning Magic kegel Master Gen2

Before Usage


Please clean Kegel Master Gen 2 with cleanser or warm water and anti-bacterial soap before first use.

After Usage

Please clean Kegel Master Gen 2 with cleanser or warm water and anti-bacterial soap and dry it.

Storage

Keep Master in a cool and ventilated place. Avoid direct sunlight and keep it away from children.

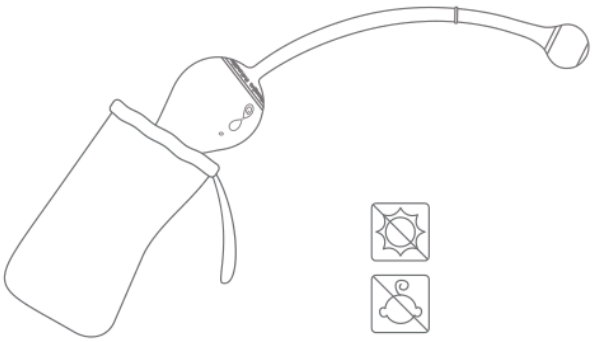Specification

Model: SMV02

Material: Body-safe silicone / ABS

Weight: 58g

Battery: 200mAh

Charging voltage: 5V

Maximum Charging Current: 500mA

Charging time: 60 minutes

Working time: 180 minutes

Waterproof: IPX6

Operating Frequency Band: 2402–2480 MHz

Maximum RF Transmit Power: −2.94 dBm

Important Safety Information

Safe Use: For adults aged 18 and over only. Use only as described in this guide. Keep away from children. Risk of injury! Remove piercings or jewellery in the genital area before use. Do not use the toy on irritated or damaged skin. Stop using the toy if you experience pain or feel unwell. Do not use the toy if it is unusually hot, damaged or deformed, or discoloured.


Charging Safety: Use the supplied charging cable and a power supply compatible with the device's rated input. Keep charging contacts dry and the cable and power supply away from water. Do not use during charging. Disconnect when fully charged.

Battery Safety: Do not disassemble, puncture or heat the device or battery, or expose them to fire. Stop use and charging if there is unusual heat, swelling, leakage or damage.

Cleaning, Storage & Disposal: Clean before and after use according to this guide, and dry before storage. Observe the stated IPX6 rating; it does not apply to the charging cable or power supply. Store dry, away from sunlight and children. Take the device and its battery to an appropriate electronics/battery collection point; do not dispose of them with household waste.


Customer Support

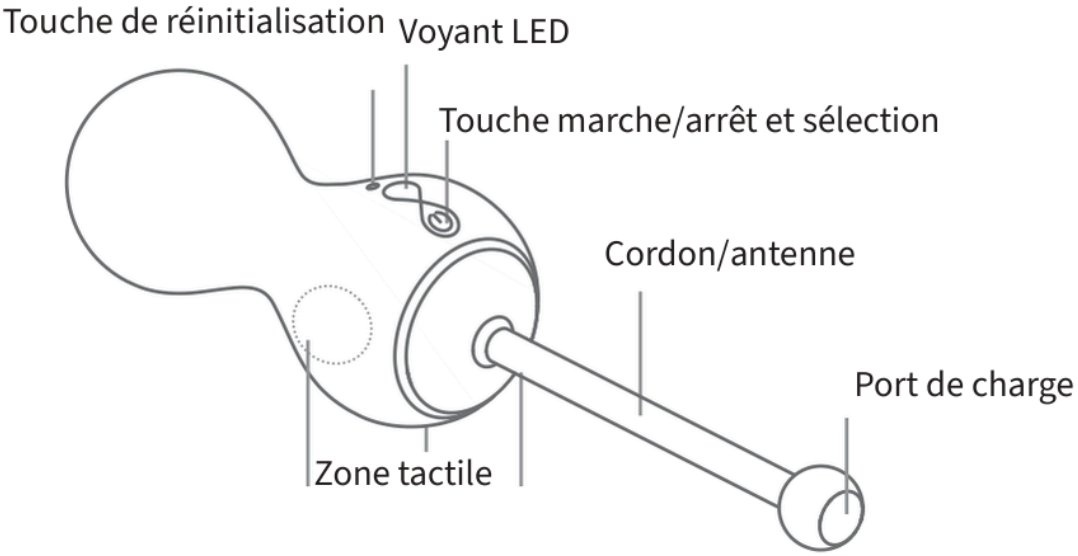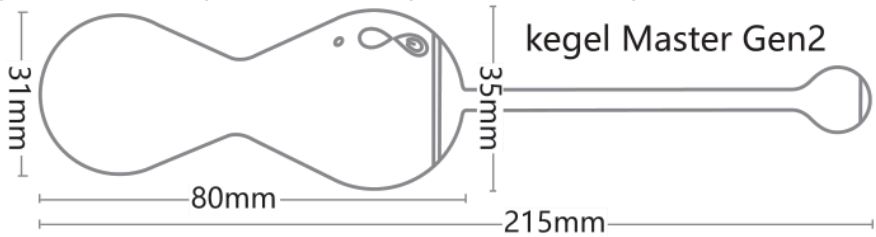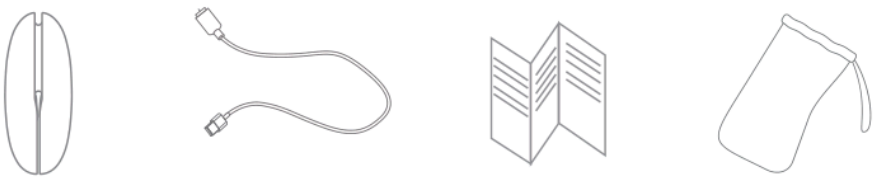
For support or warranty questions, email customercare@magicsmotion.com. Include the model and issue; include proof of purchase for warranty questions.

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Vtrump Tech (Shanghai) Co., Ltd.